





CONTINENCE CARE & UROLOGY PRODUCTS

WITH PATIENTS & SPECIALISTS WE BRING SAFETY TO UROLOGY

All together, for the future





Safe solutions for patients with chronic and acute urinary disorders

Urinary disorders affect millions of people all over the World. There are no real age limitations, no specific genders. It affects children, teenagers, adults, seniors, men, women, disabled people, and otherwise healthy people that might be treated in hospital facilities or at home.

Continence care products are meant to deal with chronic disorders for home care patients looking for a better quality of life on a daily basis for an active and independent life. Urology products are a reliable and Safe partners for health care professionals to deal with acute disorders in hospitals. For adults or chronic care patients, B. Braun provides a complete solution for safety products in terms of protocol and materials in compliance with urinary management international guidelines for health care and patient safety.




For a better understanding...


This catalogue has been designed to help identify our users.

They could be **home care patients** represented by the "big family" members according to the product, or the **health care professionals**, represented by the smaller characters here under.

In addition, a distinction between Home care and a Hospital solutions is made, choosing the small house logo for Home care & "H" logo for Hospital products.

We hope you will find our brochure instructive and enjoyable to read...

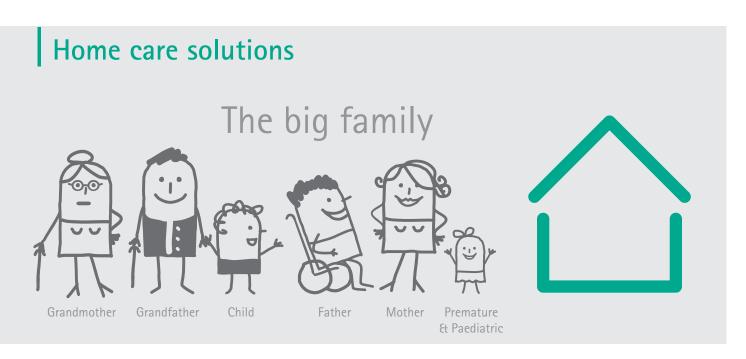

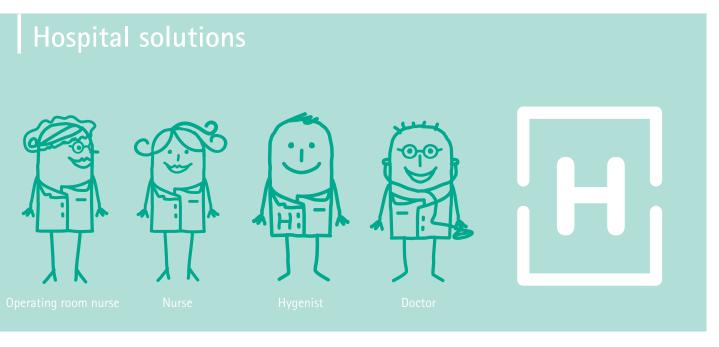

Summary



URINARY CATHETERIZATION

<u>04</u>



URINARY
COLLECTION
& MEASUREMENT

37



BOWEL MANAGEMENT & URINARY INCONTINENCE

43



WINNING COMBINATIONS

48



ORDERING INFORMATION

50





URINARY CATHETERIZATION

Intermittent Catheterization

8 ACTREEN TECHNOLOGY

More than 15 years of experience

Intermittent Catheters

10 ACTREEN® MINI RANGE (CATH & SET)

Small things can make a big difference!

12 ACTREEN® HI LITE RANGE (CATH & SET)

To use Anywhere, Anytime!

14 ACTREEN SAFESET®

The innovative solution for Urinary Tract Infection prevention - (UTI)

Suprapubic Catheters

16 CYSTOFIX®

Geared toward safety solutions

18 CYSTOFIX® SG

Allows a supra pubic catheter placement with minimal risk

Percutaneous Nephrostomy Set

19 NEPHROFIX®

Percutaneous - Nephrostomy sets

Ureteral Stents and Catheters

20 URECATH®

Ureteral stents and catheters

Urinary Dry Catheters

21 NELATON AND TIEMANN CATHETERS

Single use dry catheters for intermittent urinary catheterization

Indwelling Transurethral Catheters

22 URIMED® CATH

Silicone Foley Catheters

Indwelling Catheter Maintenance

24 URO-TAINER®

Improving Patient Quality of Life

ACTREEN® TECHNOLOGY

More than 15 years of experience


1

A unique technology based on hydropholic lubricant

Components: Water and glycerine based formula maintained troughout catheterization.


Patient's Benefits

- No mess, no splash: gel lubricant properties
- Ready to use product: no additional water required for preparation
- No risk of drying out during insertion: the level of lubrication stays the same during the insertion⁽¹⁾

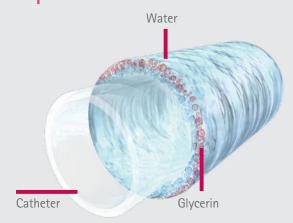



With time, most catheters dry out within the first 10min. Actreen® remains at the same level of lubrication for up to 60 minutes(*).

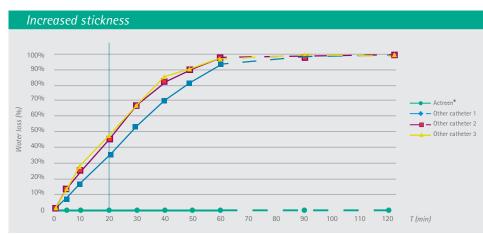

Report N° 7.5AC059.2

No risk of residue:

When the catheter is immersed in synthetic urine, Actreen® does not present any visually detectable particules in suspension⁽²⁾.

No effect on sperm mobility

ternal

An in vitro comparative study performed by an external expert demonstrates that Actreen® has negligibly toxicological effect on sperm motility and motion⁽³⁾.

⁽¹⁾ Report N°7.5AC059.2

⁽²⁾ Report N°7.5.AA.034.3

⁽a) Effect of short-term exposure to two hydrophilic-coated and one gel prelubricated urinary. Dr Auger et al. Paris, Hôpital Cochin Minerva Urologica e nefrologica, 2007;59(2):115-124

^{(4) «}Enquiry about user acceptability of Actreen® Hi-Lite Set». Data on file.

⁽⁵⁾ Report n° 7.5.AA.100.1

⁽¹⁾ For safety and hygienic reasons, we advise you to use the product in the shortest possible time after opening

The catheter features

- PVC free or DEHP free
- New material for PVC free catheter: TPO (Thermo-Plastic Polyolefin)
- Non traumatic eyes > Smooth insertion(5)







B.Braun innovative packaging No Touch technique








The catheter has no direct contact with hands and fingers

> Minimises risks of catheter contamination during handling



no touch

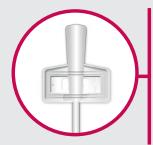

A specific connector design

> Providing an easy grip without touching the catheter(4)

Environment oriented packaging

Packaging of Actreen® is made from recyclable materials.



ACTREEN® MINI RANGE

Small things can make a big difference!

Mini catheters and sets for female urinary catheterization

Common benefits...

The Mini catheters are designed with the same concept as the other actreen products.

- Pre-lubricated with Hydrophilic lubricant, ready to use catheter, no water needed
- Safe catheter material, PVC free, DEHP free
- Non traumatic eyes(2)
- Straight forward use
- Easy to carry thanks to their low profile and lightweight





The optimal discretion for woman

The Actreen® Mini range is your perfect indoor and outdoor companion, thanks to its feminine low profile packaging and famous discreet pouch that comes with each box.







For who? For what?

Used by female patients for intermittent self urinary catheterization.

For more information

www.aotreen.com



Scan to watch instructional videos.



Brochure

⁽¹⁾ Report n°7.5.AC.067.2

⁽²⁾ Report n° 7.5.AC.020.1

⁽³⁾ Report n° 7.5.AB.043.2 (4) Report n° 7.5.AB.042.2



Actreen® Mini Cath to stay






Easy to open

- Clearly identifiable
- Easy open packaging⁽¹⁾
- Direct access to the catheter



All in one lightweight product

 Catheter + bag sealed together and folded in a low profile shape

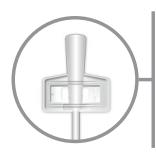

Easy to use

- Universal connector: easy to grip and handle during use
- Hygienic function: the catheter's size avoids contact with the toilets



Anti-reflux valve

 Integrated in the bag to avoid leakage during storage⁽³⁾



A feminine design


- Small and discreet packaging
- A feminine design for more discretion



Graduated bag of 700 mL

- Approximate measurement of urine
- Easy to use and empty thanks to its tear off corner⁽⁴⁾

Product de	Product description		Ø (mm)	Units/box	Reference				
ACTREEN® MINI CATH - 9 cm (catheter)									
1000		1 0	3,3	30 + 1 discreet pouch	228010*				
1		☐ 12 4,0 30 + 1 discreet pouch			228012*				
	x W	1 4	4,7	30 + 1 discreet pouch	228014*				
ACTREEN® N	MINI SET - 9 cm ((catheter)							
1000000		1 0	3,3	30 + 1 discreet pouch	239010*				
1	American Property of	□ 12	4,0	30 + 1 discreet pouch	239012*				
 9	2500	1 4	4,7	30 + 1 discreet pouch	239014*				

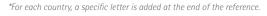



ACTREEN® HI-LITE RANGE




To use Anywhere, Anytime!

Catheter and Sets for intermittent urinary catheterization

Common benefits...

The Hi Lite range is designed with the common features below:

- Safe catheter material, PVC free, DEHP free
- Non traumatic eyes for smooth bladder drainage(1)
- Pre-lubricated with hydrophilic lubricant, ready to use catheter, no need for water
- Opening design for good and reduced dexterity
- Safe handling with two «No touch» options
- Foldable product for more discretion**

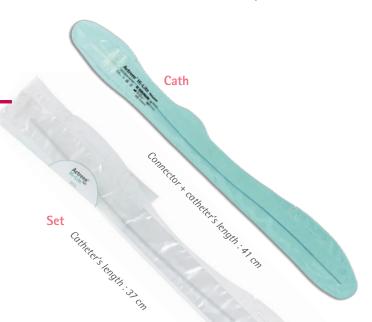

2 options to avoid touching the catheter during use.











Movable no touch

Accordion no touch

The invisible solution all day long

Its foldable shape allows to carry the product anywhere in the discreet pouch provided with each box**.









For who? For what?

Used by male, female and pediatric patients for intermittent self urinary catheterization.

For more information











Actreen® Hi-Lite cath Hi-Lite set

^{**} The CH18 and the 20 cm version can not be folded and no discreet pouch is included inside these boxes..

⁽¹⁾ Report n° 7.5.AA.100.1

⁽²⁾ Report n°7.5.AB.054.1

⁽³⁾ CE MARK in progress - Product available soon



Actreen® Hi-Lite Cath to stay





Designed to be open by everyone

- 2 options for opening adapted to the level of patient's hand dexterity
- Direct access to the catheter



All in one lightweight

- Catheter + bag sealed together and folded in a compact size**
- Fits in a pocket



Unique connector shape

• Allows to easily handle the catheter during use



Anti-reflux valve

• Integrated into the bag to avoid leakages during storage⁽²⁾



Graduated bag of 1000 mL

- Approximate measurement of urine
- Easy to read
- Bag with a friendly design
- Easy to empty thanks to its tear off corner⁽²⁾

Ordering information

Actreen® Hi-Lite Cath

	CH	ø (mm)	Units/box	Ref.
TIEMANN - 41 cm				
	8 0	2,7	30 + 1 discreet pouch	238108*
	1 0	3,3	30 + 1 discreet pouch	238110*
7	□ 12	4,0	30 + 1 discreet pouch	238112*
MAN TO SERVICE STATE OF THE SE	1 4	4,7	30 + 1 discreet pouch	238114*
	1 6	5,3	30 + 1 discreet pouch	238116*
NELATON - 41 cm				
	80	2,7	30 + 1 discreet pouch	238208*
	1 0	3,3	30 + 1 discreet pouch	238210*
P	□ 12	4,0	30 + 1 discreet pouch	238212*
miles	1 4	4,7	30 + 1 discreet pouch	238214*
	1 6	5,3	30 + 1 discreet pouch	238216*
NELATON - 20 cm				
	0 6	2,0	30 + 1 discreet pouch	238306*
	80	2,7	30 + 1 discreet pouch	238308*
	1 0	3,3	30 + 1 discreet pouch	238310*
1	□ 12	4,0	30 + 1 discreet pouch	238312*
	1 4	4,7	30 + 1 discreet pouch	238314*
l l	1 6	5,3	30 + 1 discreet pouch	238316*

Actreen® Hi-Lite Set

	CH	ø (mm)	Units/box	Ref.
TIEMANN - 37 cm				
	80	2,7	30 + 1 discreet pouch	242108*(3)
	1 0	3,3	30 + 1 discreet pouch	242110*
	□ 12	4,0	30 + 1 discreet pouch	242112*
1	1 4	4,7	30 + 1 discreet pouch	242114*
	1 6	5,3	30 + 1 discreet pouch	242116*
	1 8	6,0	30	242118*
NELATON - 37 cm				
	80	2,7	30 + 1 discreet pouch	242208*(3)
	1 0	3,3	30 + 1 discreet pouch	242210*
414	□ 12	4,0	30 + 1 discreet pouch	242212*
A	1 4	4,7	30 + 1 discreet pouch	242214*
	1 6	5,3	30 + 1 discreet pouch	242216*
	1 8	6,0	30	242218*

^{*} For each country, a specific letter is added at the end of the reference radical.

⁽³⁾ CE MARK in progress - Product available soon

ACTREEN SAFESET®

The innovative solution for Urinary Tract Infection prevention - (UTI)

Closed system concept⁽¹⁾⁽²⁾⁽³⁾ – all in one with integrated needle free sample port for single catheterization at hospital

Benefits for:

- Patient:
 - Slipperiness maintained during catheter insertion⁽⁴⁾
 - Prevents risk of contracting UTI(1,2,3)
- Personnel care:
 - Reduced risk of injury thanks to the needle free urine sample port
 - Time saving thanks to the all in one concept
- Hygenists:
 - Reliable partner for UTI reduction(1,2,3)
- Pharmacist/Purchaser:
 - Number of rationalized products and orders
 - Simplified urinary catheterization procedure and reduced waste (1)



Focus on main handling steps




Product Opening: guaranteed catheter sterility thanks to the «No touch» technique.



Desinfect the sample port before taking urine sampling,

Sterility is guaranteed inside the pouch.



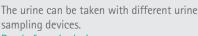

sampling devices. Ready for urinalysis.

For who?

Used by nurses in hospitals for single intermittent urinary catheterisation.

Winning



Examination and protective gloves. Vasco® Nitril long gloves.

For more information



Scan to watch instructional videos.





Brochure

Poster

- (2) ANAES (National Agency for health evaluation & accreditation) « Qualité de la pose et de la surveillance des sondes urinaires, décembre 1999»
- (3) SFHH Recommandations (French Hospital Hygienists society) «Surveiller et prévenir les infections associées aux soins, septembre 2010»
- (4) Report N°7.5AC059.2

⁽¹⁾ Infections urinaires nosocomiales de l'adulte. Conférence de Consensus SPILF & AFU, Novembre 2002, page 11

⁽⁵⁾ Report n°7.5.AB.045.1



Urinary catheterization



All in one system

 Catheter + bag integrating an urine sample port, all sealed together: no possible detachment⁽⁵⁾





Graduated bag of 1500 mL

- Approximate measurement of urine
- Tear off to empty the bag



Ready to use

 Pre-lubricated catheter, no need of extra lubricant



Anti-reflux valve

- To avoid urine backflow from the bag⁽⁵⁾
- Prevents the urine splash that could contaminate the environment



Guaranteed sterility

- «No touch» system to reduce the risk of catheter contamination
- Sterile double packaging for safe handling in operating room



Integrated needle free sample port


- Suitable with all luer sampling devices
- Designed to maintain a closed system until urine sampling⁽⁵⁾



	СН	Ø (mm)	Units/box	Reference						
ACTREEN SAFESET® NETATON – 47 cm										
	1 0	3,3	50	234210						
	12	4,0	50	234212						
	1 4	4,7	50	234214						
	1 6	5,3	50	234216						
	1 8	6.0	50	234218						

Needle free sample port

Suitable with different luer devices.

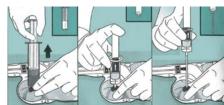

Syringe – Vaccum systems





Geared toward safety solutions

Suprapubic puncture & exchange sets

User benefits

Puncture sets:

- **Ready to use** setting: pre assembled catheter placed into the cannula, to reduce risks of sharp injury and save time
- Thin and sharp stainless steel cannula
 - 3-wings handling grip for easy prehension
- Wide range of catheter designs and sizes
 - Silicone catheters with integrated balloon for an easy replacement
 - Polyurethane pigtail catheters
 - Open tip catheters to facilitate replacement with a guide wire

Replacement sets:

- Integrated balloon silicone catheters to facilitate replacement
- Wide choice of sizes and sets to adapt to the patient's situation



Puncture sets: how to use



The catheter is already in the cannula.

Keep the protection tube on the cannula.

Pull the catheter until the tip is right below the cannula bevel.



Remove the protection tube to begin the puncture.

Winning combinations



To be used with **Urimed**® **SP** or **Urimed**® **SP**+

Replacement sets and dilation sets



Integrated balloon for easy replacement.



Short tip for patient's comfort.



Range of catheter sizes.



Range of dilators' size.



Ordering information



PUNCTION SETS WITH	I Balloon for all USE. New Range (CYSTOFIX® READY-TO-USE
--------------------	------------------------------------	------------------------

Catheter design	Catheter material	Cannula length	Cannula diameter	СН	Catheter length	Units/box	Reference
Integral balloon straight tip	Silicone	12 cm	5,4 mm	12	44 cm	10	4450200
Ready to Use	Silicone	12 cm	5,4 mm	1 4	44 cm	10	4450220
Balloon & J tip (classical balloon) Ready to Use	Silicone	12 cm	5,4 mm	12	53 cm	10	4450210

PUNCTION SETS WITH PIGTAIL CATHETER FOR SHORT TERM USE. 7 EYES. NEW RANGE CYSTOFIX® READY-TO-USE

Catheter design	Catheter material	Cannula length	Cannula diameter	СН	Catheter length	Units/box	Reference
Complete set with 2 L		8 cm	3,4 mm	1 0	65 cm	10	4450150
· ·	Polyurethane	12 cm	3,4 mm	1 0	65 cm	10	4450160
urine bag Ready to Use		12 cm	5,4 mm	1 5	65 cm	10	4450170
Foo set (ne hea)	Polyurethane	8 cm	3,4 mm	1 0	65 cm	10	4450100
Eco set (no bag) Ready to Use		12 cm	3,4 mm	1 0	65 cm	10	4450120
		12 cm	5,4 mm	1 5	65 cm	10	4450130

PUNCTION SET WITH PIGTAIL CATHETER FOR PEDIATRIC AND NEONATES
Each set contains: one split sheath cannula, protection tube, catheter. Paediatric set: + urine bag & sliding clamp. Minipaed set: + soft mandril & connectors

Catheter design	Catheter material	Cannula length	Cannula diameter	СН	Catheter length	Units/box	Reference
Paediatric complete set with 2 L urine bag Ready to Use	Polyurethane	5 cm	3,4 mm	■10	50 cm	10	4450180
Minipaed set	Polyurethane	5 cm	2,2 mm	5	50 cm	10	4441036

Catheter design	Catheter material	Cannula length	Cannula diameter	СН	Catheter length	Units/box	Reference
TUR set	Polyurethane	12 cm	5,4 mm	1 5	65 cm	10	4441354


Exchange sets and dilation sets



EXCHANGE SETS WITH INTEGRAL BALLOON CATHETER										
Set design	Catheter design	Catheter material	СН	Catheter length	Units/box	Reference				
Dilation and replacement	Balloon	Silicone	12			4450512				
with quidewire and dilator	straight tip		1 4	44 cm	5	4450514				
with guidewire and dilator	Straight tip		1 6			4450516				
		Silicone	1 0		5	4450010				
Replacement with guidewire	Balloon		12	44 cm		4450012				
Replacement with guidewire	straight tip		1 4			4450014				
			1 6			4450016				
		Silicone	12	44 cm	5	4450712				
Replacement without	Balloon		1 4			4450714				
quidewire	straight tip		1 6			4450716				
guidewire	Straight tip		1 8			4450718				
			20			4450720				
			12			4450812				
	Balloon		1 4			4450814				
Replacement single pack	straight tip	Silicone	1 6	44 cm	1	4450816				
	Straight tip		1 8			4450818				
			2 0			4450820				





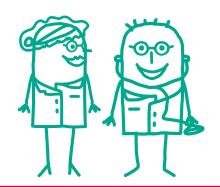

CYSTOFIX® SG

Allows a suprapubic catheter placement with minimal risk

Suprapubic sets according to the seldinger technique

User benefits

- Ready to use under local anesthesia
- Minimal pressure required to insert the 18G puncture needle in comparison with a trocar
- **Gentle dilation** of the tissues thanks to the use of a needle followed by a dilator
- Withdrawal and replacement of the catheter facilitated by the integrated silicone balloon
- Seldinger technique is well known for the insertion of Central Venous Catheter and Percutaneous Nephrostomy Catheters





18G 12 cm needle

Local anesthesia facilitated



Dilator with split sheath

• Designed for a gentle dilation of the tissues



Silicone catheter with integrated balloon

• Facilitated withdrawal and replacement thanks to the flat balloon

PUNCTION SETS WITH INTEGRAL BALLOON SILICONE CATHETERS FOR ALL USE. NEW RANGE CYSTOFIX® SG Each set contains: one scalpel, 18G needle, guidewire, dilator, catheter, 2 syringes, spigot.									
Catheter design	Catheter material	Dilator length	Dilator diameter	СН	Catheter length	Units/box	Reference		
	Silicone	15 cm	5.1 mm	1 0	44 cm	5	4450410		
Integrated balloon		15 cm	5.9 mm	12	44 cm	5	4450412		
straight tip		15 cm	6.6 mm	1 4	44 cm	5	4450414		
		15 cm	6.6 mm	1 6	44 cm	5	4450416		

Nephrostomy sets



User benefits

- Complete sets with dilators and working sheath for short and long term
- 22 cm long needle
- 1, 2 or 3 dilators for a gentle dilation of the tissues
- · Choice of pigtail and balloon catheters



Pigtail tip catheters

Ordering information

- Certon® (polyurethane)
- Low risk of incrustation
- Large drainage eyes
- Radiopaque and flexible
- Choice of 3 sizes

Balloon catheter

- Silicone
- Short tip and short balloon
- 2 eyes above balloon and 1 eye below balloon

Complete sets

	GUIDE					PIGTAIL CA	THETER		BALLOON CATHETER	
Product	Style	Tip	Ø (inch)	Ø (mm)	Length (cm)	CH8	CH11	CH14	CH12	Units/box
Nephrofix®	Seldinger	Straight	0,035"	0,9	80	4435087	4435117	4435141		5
Nephrofix® Certo	Schüller	J3	0,032"	0,8	90	4435583	4435613	4435648	4435745	5

All sets contain: one catheter, 3 pieces needle, guide wire, dilators, scalpel, syringe, ruler.

Exchange sets

	GUIDE					PIGTAIL CATHETE	R	BALLOON CATHETER	
Product	Style	Tip	Ø (inch)	Ø (mm)	Length (cm)	CH8	CH11	CH12	Units/box
Nephrofix®	Seldinger	Straight	0,035"	0,9	80	4436083	4436113	4436741	5

All exchange sets contain: one catheter, guide wire & dilator of the appropriate size.

Puncture needles: • Nephrofix®: 3 pieces 17G x 220 mm • Nephrofix® Certo: 3 pieces 18G x 220 mm - Schüller guide wire: rigid part 300 mm




URECATH®

Ureteral stents and catheters

User benefits

- Made of polyurethane (Certon®)
- X-ray visible
- Stents with smooth pigtail on bladder side designed to improve comfort

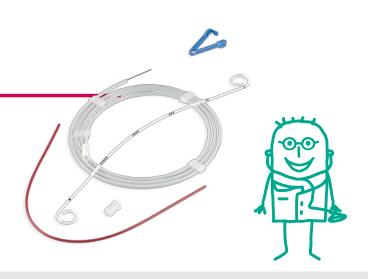

Ordering information

Stents

Stent Tip	СН		Leng	th (mm)			
Stent rip	СП	200	250	275	300	*	
OPEN / CLOSE STENT DESIGN							
	5	4430921	4430930	4430948	4430956	(1)	
Tiemann	6	4431057	4431073	4431081	4431090	(2)	
Hemann	7,5	4431103	4431111	4431120	4431138	(2)	
	9	-	4431316	-	-	(2)	
Round	7,5	4431200	4431219	4431227	-	(2)	

CH	Length (mm)						
CH	250	275	*				
OPEN /	OPEN / OPEN STENT DESIGN						
6	4437497	4437500	(3)				
7,5	4437535	4437543	(3)				
	'						

CH	Length (mm)						
СП	200	250	275	*			
CLOSE / CLOSE STENT DESIGN							
7,5	4431600	4431618	4431626	(3)			

Catheters

Tinchono	СН	Longth (om)	Number of over		Eyes location	ı	Reference	*
Tip shape	СП	Length (cm)	Number of eyes	Tip	Loop	Shaft	Reference	
URETERIC CATHE	TERS							
	4	100	4	Х			4432690	(4)
Straight	6	100	4	Х			4432703	(4)
	8	100	4	Х			4432711	(4)
	6	70	8		Х	Х	4432002	(5)
Loop	7,5	70	8		Х	Х	4432010	(5)
	9	70	8		Х	Х	4432029	(5)
Loop	7	70	4		Х		4432800	(5)

Packaging: Boxes of 10.

⁽¹⁾ Delivered with guide wire 900 mm x 0,6 mm (0,024"), pusher, cap, clamp and patient's card.

 $^{^{(2)}}$ Delivered with guide wire 900 mm x 1,1 mm (0,043"), pusher, cap, clamp and patient's card.

⁽³⁾ Delivered with guide wire 1500 mm x 0,9 mm (0,035"), 1 side curled, pusher, cap, clamp and patient's card.

⁽⁴⁾ Delivered with introducer.

 $[\]ensuremath{^{\text{(5)}}}$ Delivered with introducer, cap and clamp.



UROLOGICAL CATHETERS

Single use dry catheters for intermittent urinary catheterization

User benefits

- Medical PVC material Latex Free
- Transparent PVC allowing urine visibility throught it
- Wide range composed of Nelaton and Tiemann tips and 2 sizes for male (40 cm) and female (18 cm)
- Individually packaged sterile
- Colour coded Charrière enabling immediate size identification







The universal colour-coded connectors provide an immediate identification of the sizes.

	Length (cm)	CH	Ø (mm)	Units/box	Reference
TIEMANN CATHETER					
		■ 08	2,7	100	4500200
		1 0	3,3	100	4500201
/7		12	4,0	100	4500202
ω	40	1 4	4,7	100	4500203
l l		1 6	5,3	100	4500204
1		1 8	6,0	100	4500205
11		20	6,6	100	4500206
NELATON CATHETER					
		■ 08	2,7	100	4500207
_		1 0	3,3	100	4500208
Ph .		12	4,0	100	4500209
	40	1 4	4,7	100	4500210
1		1 6	5,3	100	4500211
		1 8	6,0	100	4500212
		20	6,6	100	4500213
NELATON CATHETER					
		08	2,7	100	4500214
P		1 0	3,3	100	4500215
	10	12	4,0	100	4500216
	18	1 4	4,7	100	4500217
T I		1 6	5,3	100	4500218
		1 8	6,0	100	4500219



URIMED® CATH

Silicone Foley Catheters

Indwelling balloon catheters

User benefits

Silicone catheters are hypoallergenic⁽¹⁾ and therefore preferable for long-term settings:

- Less tissue irritation
- Reduced risk of encrustation in long-term catheterized patients⁽¹⁾
- Less traumatic catheterization and less risk of blockage compared to latex thanks to thinner walls and larger lumen, for the same charrière size(1)(2)
- Visible on X-Rays




Double sterile packaging
2 openings: lateral and longitudinal



Wrinkle the primary peel-pack.



Tear it off with the side pre-cuts...

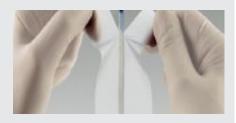

or with the pre-cuts down the catheter.

Winning combination

 Urine Measurement Ureofix® 500 Classic & Compact



- Urine Collection 2 Way: Urimed® SP bags
- 3 Way Prostatic catheters: 4 L bag



 Bladder washout, irrigation and instillation Uro-Tainer®



For more information



Brochure



Adapted to every need



2 Way Nelaton

- 1st intention product
- Suitable for both Male & Female
- For routine bladder drainage
- Closed tip
- 2 drainage eyes



2 Way Tiemann


- Designed to ease catheterization in Male
- For routine bladder drainage
- Closed tip
- 1 drainage eye



3 Way Dufour Prostatic Catheters

- A 3rd channel enables irrigation after prostate surgery
- For routine bladder drainage and irrigation and/or instillation
- Open tip
- 3 drainage eyes + 1 bevel

	CH	Balloon capacity (mL)	Units/box	Reference
URIMED® CATH FOLEY N	ELATON			
	12	10	10	4563250
	1 4	10	10	4563251
	1 6	10	10	4563252
	1 8	10	10	4563253
	2 0	10	10	4563254
~	2 2	10	10	4563255
	2 4	10	10	4563256
URIMED® CATH FOLEY TI	EMANN			
	12	10	10	4563257
	1 4	10	10	4563258
	1 6	10	10	4563259
and the same of th	1 8	10	10	4563260
~	2 0	10	10	4563261
URIMED® CATH FOLEY D	UFOUR			
	1 6	40	10	4563262
**************************************	1 8	50	10	4563263
	2 0	60	10	4563264
delle	2 2	70	10	4563265
	2 4	80	10	4563266

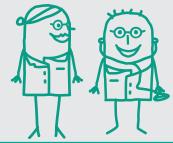

⁽¹⁾ Geng V et al. Evidence-based Guidelines for Best Practice in Urological Health Care: Catheterisation - Indwelling catheters in adults - urethral and suprapubic. EAUN, 2012

⁽²⁾ Morris NS et al. Which indwelling urethral catheters resist encrustation by Proteus mirabilis biofilms? BJU 1997; 80:58-63

URO-TAINER®

Improving Patient Quality of Life

User benefits

Uro-Tainer® technology

- Patients with indwelling urinary catheters are at high risk of developing complications such as infection and encrustation
- There is a strong rationale for on-going catheter maintenance to prevent complications and maximize the life of the catheter together with patient comfort and quality of life
- For over 25 years the Uro-Tainer® technology and solutions have been used in clinical practice to manage encrustation and debris formation

Uro-Tainer® solutions

• TWIN:

Prevention & removal of blockage

• Uro-Tainer® Polihexanide:

Bacterial decolonisation

• Uro-Tainer® NaCl:

Bladder instillations



For who?

50%

of catheterized patients affected due to encrustation and blockage of indwelling urinary catheters⁽¹⁾.

Winning combination

Urimed® Cath Silicone Foley Catheters



For more information



Brochure

^[1] Getliffe KA. The characteristics and management of patients with recurrent blockage of long-term urinary catheters. J Adv Nurs 1994; 20: 140-149



Uro-Tainer® Solutions





Prevention & removal of blockage



Bacterial decolonisation



Bladder instillations

Twin

• Suby G

Indication: A mildly hypotonic fluid that is less irritatinng as a result of the addition of magnesium. This fluid is specially designed to prevent phosphate crystallization and dissolve existing calcification in indwelling catheters.

• Solutio R

Indication: A mildly hypotonic fluid that is specifically designed for cathethers with stubborn cacification where Suby G does not provide a sufficient result. This is due to its highter concentration of citric acid and gluconolactone.

In addition, this fluid minimises trauma when removing an indwelling catheter.

Uro-Tainer® Polihexanide

Indication: Uro-Tainer® Polihexanide 0.02% is used for routine decolonisation (bacteria removal) of the catheter.

Uro-Tainer® NaCl

• NaCl o.9%

Indication: This isotonic fluid is used primarily for mechanically cleaning the bladder and catheter.
eg. in the case of debris forming in the bladder.

• M NaCl 0.9% with injection port Indication: The Uro-Tainer® M has an injection port that is specifically for administering drugs.

Products manufactured by: B.Braun Medical AG Seesatz 17 6204 Sempach Switzerland

Ordering information

URO-TAINER® Twin

	Size	Units/box	Reference					
URO-TAINER® TWIN SUBY G								
Standard	100 mL	10	FB99839					
Twin	60 mL	10	9746609					
URO-TAINER® TW	URO-TAINER® TWIN SOLUTION R							
Standard	100mL	10	FB99841					
Twin	60 mL	10	9746625					

URO-TAINER® Polihexanide

	Size	Units/box	Reference			
URO-TAINER® POLIHEXANIDE						
	100 mL	10	FB99965			

URO-TAINER® NaCl

	Size	Units/box	Reference				
URO-TAINER® NACL 0.9%							
	100 mL	10	FB99833				
URO-TAINER® M NACL 0.9% with injection port							
	100 mL	10	FB99853				







URINARY COLLECTION & MEASUREMENT

Urine Measurement Systems

28 UREOFIX® 500 CLASSIC

Accurate and safe urine output management

30 UREOFIX® COMPACT

Accurate and safe urine output management

Closed Urine Drainage Bags

32 UREOFIX® 112 PLUS

Closed urine drainage system with triple safety concept

33 POST OPERATIVE 4 L BAG

Closed urine drainage system with large capacity for post TUR use

34 URIMED® SP

Quality within everyone's reach

Adult Urine Bags

36 URIMED® BAGS

A sophisticated range of urine bags for daily use

37 URIMED® ACCESSORIES

A complementary range of accessories

38 URIMED® B'BAGS

Adult urine collection bag

Pediatric Urine Bags

40 URINOCOL®

Looking for adapted solutions

41 PEDIATRIC URINE BAGS®

Cost effective range for urine pediatric collection

UREOFIX® 500 CLASSIC

Accurate and safe urine output management

3 chambers Urine Measurement System

User benefits

- Precise measurement with 3 chambers
- For reduced risk of retrograde migration of bacteria thanks to the Pasteur chamber and 2 non return valves
- Safety for the user thanks to the needle free sample port
- **User convenience** with a wide choice of bags sizes and tube lengths
- **Direct** transfer of urine from chamber to bag with the drainage tab





User convenience



Needle Free sample port



Clear marking drainage tab



Cross drainage outlet with dog house



Technical aspects



Closed System

- Dry Pasteur chamber
- With non return valve



Non return valve

• Between chamber and bag



Accurate measurement

- Precise measuring chamber up to 50 mL
- 3 chambers for precision in all situation



Choice of bag

- 2 L bottom outlet bag
- 1.5 L & 3.5 L closed bag



Choice of fixation

- 2 straps for permanent fixation
- 1 built-in hook
- 1 pin for use of optional pin strap



Winning combination

• Urimed® Cath Silicone Foley Catheters

	lubing length (cm)	Units/box	Reference
UREOFIX® 500 CLASSIC SETS WITH CLOSED BAGS			
Ureofix® 500 Classic with exchange bag 1,5 L	120	10	4417910
Ureonx 500 classic with exchange day 1,5 L	170	10	4417920
Ureofix® 500 Classic with exchange bag 3,5 L	170	10	4417960
UREOFIX® 500 CLASSIC SETS WITH BOTTOM OUTLET BAGS			
	120	10	4417930
Ureofix® 500 Classic with premounted b.o bag 2 L	170	10	4417940
	200	10	4417950
UREOFIX® 500 CLASSIC SPARE BAGS	Volume		
Livertiv® FOO Classic classed have	1.5 L	100	4417535
Ureofix® 500 Classic closed bag	3.5 L	100	4417543
Ureofix® 500 Classic Bottom Outlet Bag	2 L	100	4417551
UREOFIX® 500 CLASSIC ACCESSORIES			
Fixation strap		10x10	4417755



UREOFIX® COMPACT

Accurate and safe urine output management

2 chambers compact Urine Measurement System

User benefits

- Precise measurement with 2 separated chambers
- For reduced risk of urinary tract infections thanks to the Pasteur chamber and 2 non return valves
- Safety of use thanks to the needle free sample port
- Convenience of use with 3 practical fixation means and the compact size
- Quick transfer of urine from chamber to bag. (5" for 80 mL - Test report 7.10W019.1)





User convenience



Double hooks for quick and stable fixation



2 straps for permanent fixation



Pin strap for permanent fixation and quick repositioning



Needle free sample port to avoid sharp injury



Technical aspects



• Pasteur chamber with non return valve



• Non return valve between chamber and bag



• Large diameter tube between chamber and bag for quick urine transfer




 Accurate urine measurement. Precise and easy to read measuring scale in the critical range 0-50 mL



• Compact design to facilitate the placement on the bedframe and avoid contamination by the floor



• Urimed® Cath Silicone Foley Catheters

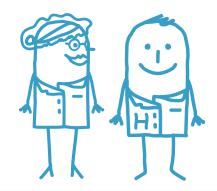

Production description	Tubing/Lenght (cm)	Units/Box	Reference
UREOFIX® COMPACT			
UREOFIX® Compact - 2 Chambers - 2 L bag	150	10	4417901
ONEOFIX® Compact - 2 Chambers - 2 L bag	120	10	4417902
Fixation strap	-	10x10	4417755





UREOFIX® 112 PLUS

Closed urine drainage system with triple safety concept


User benefits

Triple safety concept:

- Infection prevention thanks to the drip chamber and the two non return valves
- Safety for the users with the needle free sample port
- **Bottom outlet stored** in a vented protective cap to avoid contamination with the floor
- Convenient fixation means with 2 hooks and 1 cord



User convenience

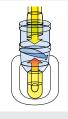







Safe drip chamber

Designed to reduce the risk of retrograde bacteria migration. Air-Vent for a pressure balance during drainage.

Double Safety System:

Anti-Reflux valves within the Drip Chamber and within the bag.

Needle free sample port

Latex-free with universal connector.

Eliminates the risk of needle stick injury during urine sampling.

Vented protective cap on the bag

For cross outlet fixation to avoid contamination.

Ordering information Units/Box Reference UREOFIX® 112 PLUS 25 284315



+ Features:

- 2000 mL bag capacity Sterile
- Graduated scale from 25 to 2000 mL
- Tube of 110 cm
- 2 hanging systems (rope hanger + universal double-hook for bed fixation)
- Cross outlet with stopper
- Latex free gloves included



POST OPERATIVE 4 L BAG

Closed urine drainage system with large capacity for post TUR use

User benefits

- 4 L capacity for use after TURP and TURB
- **Prevention of infection** thanks to the drip chamber and the 2 non return valves
- Safety for the users with the needle free sample port
- **Bottom outlet stored** in a vented protective cap to avoid contamination with the floor
- Convenient fixation means with 2 hooks and 1 cord



User convenience



Needle free sample port on the staged connector. Eliminates the risk of needle stick injury during urine sampling.



Tube without and with pump

- Standard model without pump for use with high irrigation flow
- Model with pump to allow washouts without disconnecting the catheter



2 non return valves: Anti-Reflux valves in the Drip Chamber and within the bag.



Winning combination

To be used with the 3 way silicone catheter Urimed® Cath Dufour.

Ordering information

Production description	Units/Box	Reference
4 L POST OPERATIVE URINE BAG		
4 L bag TUR*	20	4801020
4 L bag TUR with pump*	20	4801030

*Delivered with clamp, bed sheet clip in double pouch.



URIMED® SP

Quality within everyone's reach

Urine bags with needle free sample port.


User's benefits

- Reduced risk of infections thanks to the integrated anti-reflux valve
- Safety of use with the needle free sample port
- Cross contamination with the floor is limited thanks to the outlet storage on the bag
- Non kinking tube with large diameter for a better flow of urine
- Universal connection



Urimed® SP+

Focus on main handling steps

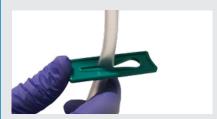

Clamp the tube in order to stop the urine flow



Desinfect the needle free sampling port located on the connector



Take the urine from the needle free sample port with a luer outlet device such as a syringe needle free

For whom? For what?

Used by nurses for urine collection and sampling during a short term period of time.

Winning combination



CYSTOFIX® Ready to use Suprapubic puncture sets.

For more information



Brochure





Technical characteristics



Clamp

• Used to stop urine flow during sampling



Needle Free Sample Port

• Eliminates the risk of needle sharp injury during urine sampling



Air filter

- Allows a complete filling of the bag
- Provides additional safety when the pressure is applied on the bag



Cross outlet

- Easy to manage, with only one hand
- Free flow, avoids urine stagnation and limits risk of infection



Outlet Storage

 The outlet can be stored to avoid contamination with the floor risk



Additional feature on URIMED SP+

Double Hanging System:

- Universal double-hook for bed fixation
- Rope hanger facilitating patient movement

	Length tube (cm)	Ø Tube (mm)	Units/box	Reference
URIMED® SP 2 L				
	100	9	10	28610
	100	9	50	28650
URIMED® SP+ 2 L				
O	130	9	30	28630



URIMED® BAG

A sophisticated range of urine bags for daily use

Adult urine collection bag

Patient benefits

- A highly sophisticated product range for a large range of application
- Designed to be easy to use
- Safe usability thanks to the anti-reflux valve
- Secure drainage: outlets designed to be easy to handle
- Non kinking tube with large diameter for a better flow of urine
- Universal connection



Main features

Night and Day bags (2 L)



Large capacity of urine collection (2 L)



Free flow cross outlet with stopper



A soft back foil on the bags to avoid skin irritation



Free flow lever tap designed to be easy and secure to drain with one hand

Urimed® Tribag Plus 800 mL





Specific 3 chamber pouch for a comfortable leg fit.

Ordering information

Night and Day bags (2 L)

	Length tube	Ø Tube (mm)	Units/box	Reference	
URIMED® BAG 2 L					
Non sterile	120 cm connected connector	9	10	28300A	
MEDICARE® 2 L					
Sterile	90 cm	6.5	10	28115	
Non sterile	90 cm	6.5	10	28114	

Leg bags

	Length tube	Ø Tube (mm)	Units/box	Reference	
URIMED® BAG PLUS 0,5 L					
Non sterile	60 cm cut to fit	9	10	28501A	
URIMED® TRIBAG PLUS 0,8 L					
3 Chamber – Sterile	60 cm connected connector	9	10	28305	
Non sterile	60 cm cut to fit	9	10	28306	
URIMED® BAG PLUS 1.5 L					
Non sterile	90 cm cut to fit	9	10	28150A	



URIMED® ACCESSORIES

A complementary range of accessories

Patient benefits

- A useful range of accessories to use in complement with the urine bags
- Practical and adapted to different applications





Urimed® Fix

 Can be worn on the thigh or lower leg



Urimed® Straps

• Adapted for leg bag



Urimed® Klett

• Fixation for tube and catheter



Urimed® Handle

To hang at the bedside or wheelchair

Ordering information

Accessories

	Tube	Units/box	Reference	
URIMED® FI	X*			
	Small	2	68520R	
	Medium	2	68530R	
	Large	2	68540R	
URIMED® STRAPS				
		1 pair	68550R	

*	Distributed	by	В.	Braur

	Tube	Units/box	Reference		
URIMED® K	LETT				
	Small	10	28260		
	Large	10	28261		
URIMED® HANDLE					
		1	68548R		

URIMED® B'BAGS

Adult urine collection bag

User benefits

- Wide range of urine bags:
 - Closed and Drainable bags
 - Large capacity: 2 L and 1,5 L
 - Sterile Non Sterile
 - With or without Anti-reflux valve
- Latex free material
- Graduated scale on all bags
- Patient data can be filled in the multilingual table
- Universal connection



Types of bags

Closed bags



For collection of urine and secretion

Drainable bags



• For urine collection

For who? For what?

Used by the nurses or by the patients for urine collection and secretion.

Winning combination



Urimed ® Cath
Indwelling balloon catheters,
Silicone Foley Catheters

For more information



Brochure





Closed bags



Tear off for emptying the bag.



Drainable bags



Free flow cross outlet designed to be managed with only one hand.



Outlet can be stored to limit risk of contamination.



Holes for universal hanging system



Anti-reflux system

Closed bags exist with and without anti-reflux valve and all drainable bags include an anti-reflux valve to avoid urine backflow during use.

Ordering information

Closed bags

		Capacity (L)	Units/box	Reference		
URIMED® B'BAG – Tube length: 90 cm – Tube diameter intR : 5.05 & extR 6.5						
Sterile	No ARV*	2	100	29400		
bags	NO ANV	1,5	100	29420		
	No ARV* With ARV*	2	250	29415		
Non sterile		1,5	250	29435		
bags		2	250	29410		
		1,5	250	29430		
Non sterile bags	No ARV*	2	30 + 1 hanging accessory**	29416		
Non sterile bags	With ARV*	2	30 + 1 hanging accessory**	29411		





**Hanging accesory

Some references exist with an hanging accessory incuded in the box.

Drainable bags

	_			
		Capacity (L)	Units/box	Reference
URIMED® B'BAG - Tube length: 90 cm - Tube diameter intR : 5.05 & extR 6.5				
Sterile bags	With ARV*	2	100	29440
Non sterile	With ARV*	2	200	29450
bags	VVILII ANV	1,5	200	29455
Non sterile bags	With ARV*	2	30 + 1 hanging accessory**	29451






URINOCOL®

Looking for adapted solutions

Pediatric urine bags

Patient benefits

- Wide sterile range designed to be adapted to all babies including premature ones
- Shape of adhesive adapted to the boy and girl anatomy
- High skin tolerance adhesive: gentle on application and removal
- Safe material: PVC Free / Latex free
- Permanent visual control of urine thanks to the transparent material
- Bags available with or without drain and anti-reflux valve

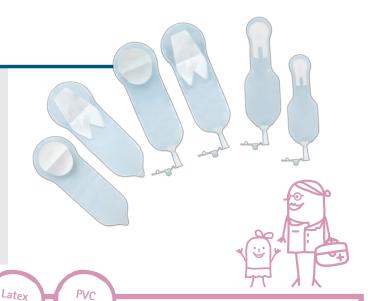

Main features

For baby > 2500g

Designed to fit the baby anatomy



One specific shape for boys, one for girls

Suitable for different indications





Sampling with syringe

Connection to urine bag





Depending on the baby weight – exists in 2 sizes: Standard et Mini

Parachute shape

Designed to maintain the bag in place more effectively.

Available with extension line: adapted for incubator.

Weight (g) Units/box Reference

1. No drainage Closed bag

Urine sampling for biochemical, cytological or bacteriological examination

2. With drainage Opened bag collecting bag for diuresis management

Ordering information

		UTITES/ UUX	nercrence			
URINOCOL PAEDIATRIC						
Sterile open	Boy	50	227546*			
paediatric bags	Girl	50	227566*			
Sterile closed	Boy	100	227550*			
paediatric bags	Girl	100	227560*			

*For each country, a specific letter is added at the end of the reference.

Units/	DOX	Reference

			5 (5)				
URINOCOL PREM	URINOCOL PREMATURE						
	Sterile	Standard	1200-2500	50	227571*		
}	premature bags	Mini	<1200	50	227570*		
	Sterile premature bags	Standard	1200-2500	50	227573*		
1	(with extension line**)	Mini	<1200	50	227572*		

**made of PVC: DEHP free



STANDARD PEDIATRIC URINE BAGS

Cost effective range for urine pediatric collection

User benefits

- Self adhesive unisex shape
- Size of the stick surface is variable
- Graduated scale at 10 mL up to 100 mL
- Bags available with or without drain and anti-reflux valve
- Transparent material film









- A unisex shape
- Self adhesive, the size of the stick surface is variable for boys and girls

 A graduated scale on the bag up to 100 mL to easily measure the urine throught the film (approximate measurement)

Ordering information - To be confirmed

	Sterile / Non sterile	ARV*/ Without ARV*	Closed / Drainable	Units/box	Reference
STANDARD PEDIAT	RIC URINE BAGS				
	Non sterile	Without ARV*	Drainable - With drainage tube of 90 cm - With slide clamp	20 – Individually packed	4423046
For Children	Non sterile	With ARV*	Drainable - With bottom outlflow nozzle	10 - Bulk packed	4423135
	Sterile – Individually wrapped	Without ARV*	Closed	10 - Individually packed	4425030
	Non sterile	Without ARV*	Closed	10 - Bulk packed	4423038
For Neonates	Sterile – Individually wrapped	Without ARV*	Closed	10 - Individually packed	4425049





BOWEL MANAGEMENT URINARY & INCONTINENCE

Colonic Lavage - Colonic Irrigation

44 RECTOBAG®

Enema and Diagnosis Closed Systems

Accessories

45 RECTAL CATHETERS

For enema purposes and intestine evacuation or for the injection of contrast media

Male External Catheters

46 URIMED® VISION

Connect to freedom

Female Intravaginal Tampon

47 DIVEEN®

Feel free again



RECTOBAG®

Enema and Diagnosis Closed Systems

User benefits

- Single use closed enema kit always at hand
- 3 different kits available
- **Hygienic handling:** Rectal balloon catheter designed to be lodged in the rectum
- 2 L bag capacity with eyelets
- 100% silicone rectal balloon catheter
- Latex-Free





 2 L lavage bag to be closed with liquid-tight cap



• Silicone soft balloon rectal catheter



Special set for X-Ray examination procedures

	Product description		Reference
1/1000	Closed System for intestinal lavage • 2 L lavage bag • Rectal tube with balloon • 150 cm large bore tubing with side clamp Sterile / single-packed	20	4423364
	Closed System for intestinal lavage and irrigation 2 L lavage bag 150 cm large bore tubing with side clamp Y-piece with connector and side clamp Sterile / single-packed	je 20	4423348
	Closed System for lavage bag • 2 L lavage bag • Connector Non sterile / bulk packed	40	4423356



RECTAL CATHETERS

For enema purposes and intestine evacuation or for the injection of contrast media

User benefits

- Latex-free
- Individually packaged in box of 100
- Nelaton 40 cm length
- Immediate size identification
- Colour coded Charrière enabling immediate size identification







 The universal colour-coded connectors provide an immediate size identification.



Length (cm)	CH	Units/box	Reference
	■ 22	100	4500300
40	■ 25	100	4500301
40	■ 28	100	4500302
	■ 30	100	4500303



URIMED® VISION

Connect to freedom

100% Silicone Male External Catheter

User benefits

Male External Catheters offer a discreet and reliable solution to male urinary incontinence.

The real benefit for patients is to continue having an active life with minimal inconvenience in a secured way.


- Security
 - Skin-friendly: excellent skin tolerance thanks to the silicone material
 - Silicone material enables high performance in water vapour permeability to avoid maceration

Autonomy

- Connection to all major types of collecting bags
- Connected to urine bag, Urimed® leg & bed bags

Convenience

- One piece system

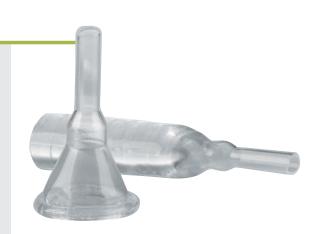



Kink proof silicone funnel

• Helps to ensure uninterrupted urine flow

Perfect Fit! Choose the right connection!

	Ø Tube (mm)	Units/box	Reference		
URIMED® VISION STANDARD - 7,5 cm length with 4,5 cm adhesive area					
	25	30	IH2525		
	29	30	IH2529		
	32	30	IH2532		
	36	30	IH2536		
	41	30	IH2541		
	25	30	IH3525		
1	29	30	IH3529		

41	30	IH2541
25	30	IH3525
29	30	IH3529
32	30	IH3532
36	30	IH3536
41	30	IH3541
	- 7,5 cm length wi 25 29 32 36	- 7,5 cm length with full adhesive 25 30 29 30 32 30 36 30

	Ø Tube (mm)	Units/box	Reference
URIMED® VISION SHORT			
	25	30	IH4525
	29	30	IH4529
	32	30	IH4532
	36	30	IH4536
	41	30	IH4541

Description	Units/box	Reference
Urimed® Bag Plus 500 mL	10	28501A
Urimed® Bag Plus 800 mL	10	28306A
Urimed® Bag Plus 1500 mL	10	28150A
Urimed® Bag 2000 mL	10	28300A







DIVEEN®

Feel free again

For stress or mixed urinary incontinence(1)

Diveen® is a intravaginal device indicated for women suffering from stress or mixed urinary incontinence.

What is Diveen®?

As simple as a tampon to apply!⁽¹⁾


B.Braun patent⁽¹⁾









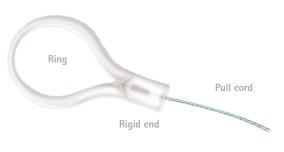



User benefits

- A unique mechanism enabling signicant reduction of the risk of urine loss.⁽¹⁾
- 100% discreet protection.
- Biocompatible material⁽¹⁾.
- 12 hours of incontinence control to enjoy physical activity with total peace of mind. (2)
- Easy to use. No need to remove Diveen® during the day, even when going to the toilet.
- Can be inserted without applicator.

Innovative and patented concept⁽¹⁾

Diveen® has been designed to support the anterior vaginal wall in the area below the bladder neck. During an effort, the weight of the bladder puts pressure on the flexible ring of the device. This causes the rigid part to move and ensures closure of the urethra, thus signicantly reducing the risk of urine loss.⁽²⁾

Clinical study

Clinically proven efficacy in the management of stress urinary incontinence in women. (2a)

• Diveen® reduces 69% of urine loss on average. (2a)

		Size	Units/box	Reference
DIVEEN®				
050	Take.	Small	5 Devices + 1 Applicator	3132
1	10 39 10	Medium	5 Devices + 1 Applicator	3133
and the same	2	Small	15 Devices + 1 Applicator	3022
	-	Medium	15 Devices + 1 Applicator	3023

⁽¹⁾ Patent EP 1 924 219 B1

⁽²⁾ Cornu J.N. et al. 75nc007 device for noninvasive stress urinary incontinence management in women: a randomized controlled trial (55 patients were controlled and analyzed - 26 controls and 29 treated). Total of 41 analyzed (24 in the control group and 17 in the treatment group). International Urogynecology Journal 2012, 23(12):1727-34.

⁽²⁰¹⁾ After 14 days, reduction of incontinence episodes frequency (primary parameter) of -68,8% ± 22% versus -8,25% ± 25,4% in the control group (p<0,0001) in per-protocol analysis.















WINNING COMBINATIONS

50	CYST0FIX®
	or
	URIMED® SP+ 2 L
	WITH BED
	FIXATION
	Vasco® Nitril long gloves - Sterile Hand protection Softaskin®
1/10	le External Catheters
50	URIMED® VISION
	URIMED® TRIBAG
	or
	URIMED® BAG
	1.5 L AND 2 L
	Urimed® Handle – Urimed® Straps
Ind	welling and Prostatic Catheters
50	URIMED® CATH
	or
	POST OPERATIVE 4 L BAG
	or
	UREOFIX® COMPACT
	or
	or UREOFIX® 112 or
	or UREOFIX® 112


CREATE YOUR SOLUTION

To make it concrete, B. Braun is proposing solutions called «The winning combinations».

«Winning» in terms of time saving, efficiency in the results, safety for patients & health care professionals.

Suprapubic Catheters

CYSTOFIX®

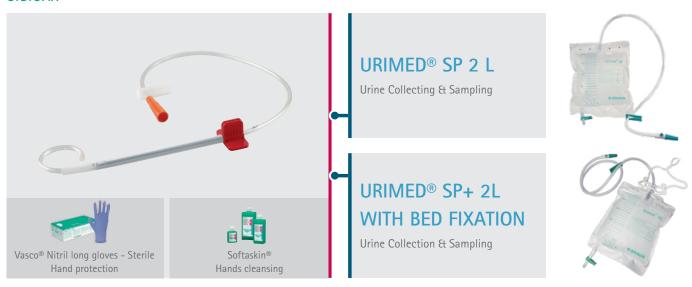

Male External Catheters

URIMED® VISION



Indwelling Transurethral Catheters

URIMED® CATH



ORDERING INFORMATION





ACTREEN® MINI RANGE

Product description	CH	Ø (mm)	Units/box	Reference
ACTREEN® MINI CATH - 9 cr	n (catheter)			
Y000	1 0	3,3	30 + 1 discreet pouch	228010*
	□ 12	4,0	30 + 1 discreet pouch	228012*
	1 4	4,7	30 + 1 discreet pouch	228014*
ACTREEN® MINI SET - 9 cm	(catheter)			
167/	1 0	3,3	30 + 1 discreet pouch	239010*
F.	□ 12	4,0	30 + 1 discreet pouch	239012*
	1 4	4,7	30 + 1 discreet pouch	239014*

^{*}For each country, a specific letter is added at the end of the reference.



ACTREEN® HI-LITE RANGE

Actreen® Hi-Lite Cath

	CH	ø (mm)	Units/box	Ref.
TIEMANN - 41 cm				
	■ 08 ■ 10 □ 12 ■ 14 ■ 16	2,7 3,3 4,0 4,7 5,3	30 + 1 discreet pouch 30 + 1 discreet pouch 30 + 1 discreet pouch 30 + 1 discreet pouch 30 + 1 discreet pouch	238108* 238110* 238112* 238114* 238116*
NELATON - 41 cm	10	5,5	30 + 1 discrect pouch	230110
	■ 08 ■ 10 □ 12 ■ 14 ■ 16	2,7 3,3 4,0 4,7 5,3	30 + 1 discreet pouch 30 + 1 discreet pouch 30 + 1 discreet pouch 30 + 1 discreet pouch 30 + 1 discreet pouch	238208* 238210* 238212* 238214* 238216*
NELATON - 20 cm				
	■ 06 ■ 08 ■ 10 □ 12 ■ 14 ■ 16	2,0 2,7 3,3 4,0 4,7 5,3	30 + 1 discreet pouch 30 + 1 discreet pouch	238306* 238308* 238310* 238312* 238314* 238316*

Actreen® Hi-Lite Set

	CH	ø (mm)	Units/box	Ref.
TIEMANN - 37 cm				
	8 0	2,7	30 + 1 discreet pouch	242108*(3)
	1 0	3,3	30 + 1 discreet pouch	242110*
	□ 12	4,0	30 + 1 discreet pouch	242112*
9	1 4	4,7	30 + 1 discreet pouch	242114*
	1 6	5,3	30 + 1 discreet pouch	242116*
	1 8	6,0	30	242118*
NELATON - 37 cm				
	8 0	2,7	30 + 1 discreet pouch	242208*(3)
	1 0	3,3	30 + 1 discreet pouch	242210*
1	□ 12	4,0	30 + 1 discreet pouch	242212*
A	1 4	4,7	30 + 1 discreet pouch	242214*
	1 6	5,3	30 + 1 discreet pouch	242216*
	1 8	6,0	30	242218*

^{*} For each country, a specific letter is added at the end of the reference radical.

[3] CE MARK in progress - Product available soon

H AC

ACTREEN SAFESET®

	CH	Ø (mm)	Units/box	Reference
ACTREEN SAFESET® NETATO)N - 47 cn			
	1 0	3,3	50	234210
	□ 12	4,0	50	234212
	1 4	4,7	50	234214
4	1 6	5,3	50	234216
	1 8	6,0	50	234218

URINARY CATHETERIZATION

Suprapubic Catheters



CYSTOFIX®

п.		_4.			L_
М	ın	CTI	ıre	set	rs.
	•••	~~~		30	-



PUNCTION SETS WITH BALLO Each set contains: one split sheat					oper.		
Catheter design	Catheter material	Cannula length	Cannula diameter	СН	Catheter length	Units/box	Reference
Integral balloon straight tip	Silicone	12 cm	5,4 mm	12	44 cm	10	4450200
Ready to Use	Silicone	12 cm	5,4 mm	1 4	44 cm	10	4450220
Balloon & J tip (classical balloon) Ready to Use	Silicone	12 cm	5,4 mm	12	53 cm	10	4450210

Catheter design	Catheter material	Cannula length	Cannula diameter	CH	Catheter length	Units/box	Reference
Complete set with 2 L		8 cm	3,4 mm	10	65 cm	10	4450150
urine bag	Polyurethane	12 cm	3,4 mm	1 0	65 cm	10	4450160
Ready to Use		12 cm	5,4 mm	1 5	65 cm	10	4450170
Eco set (no bag)		8 cm	3,4 mm	1 0	65 cm	10	4450100
` 5'	Polyurethane	12 cm	3,4 mm	1 0	65 cm	10	4450120
Ready to Use		12 cm	5.4 mm	15	65 cm	10	4450130
UNCTION SET WITH PIGTAIL (ATHETER FOR PENIATRIC		5,4 11111	13	05 CIII	10	7730130
PUNCTION SET WITH PIGTAIL (ach set contains: split sheath can		AND NEONATES				10	++30130
		AND NEONATES				10	4450180
ach set contains: split sheath can aediatric complete set with 2 L urine bag	nula, protection tube, cathet	CAND NEONATES er. Paediatric set: + urine l	oag & sliding clamp. Minipaed	set: + soft m	andril & connectors.		

Exchange sets and dilation sets



EXCHANGE SETS WITH INTEGRAL BALLOON CATHETER								
Set design	Catheter design	Catheter material	CH	Catheter length	Units/box	Reference		
Dilation and replacement			12			4450512		
	Balloon straight tip	Silicone	1 4	44 cm	5	4450514		
with guidewire and dilator			1 6			4450516		
			1 0			4450010 4450012 4450014 4450016 4450712 4450714		
Replacement with guidewire	Balloon straight tip	Silicone	12	44 cm	5	4450012		
	Danoon Straight tip	JIIICOIIC	1 4	44 (111	3	4450014		
			1 6			4450016		
			12			4450712		
Replacement without			1 4			4450714		
quidewire	Balloon straight tip	Silicone	1 6	44 cm	5	4450716		
guidewire			1 8			4450718		
			2 0			4450720		
			12			4450812		
			1 4			4450814		
Replacement single pack	Balloon straight tip	Silicone	1 6	44 cm	1	4450816		
			1 8			4450818		
			2 0			4450820		

Н CYSTOFIX® SG

PUNCTION SETS WITH INTEGRAL BALLOON SILICONE CATHETERS FOR ALL USE. NEW RANGE CYSTOFIX® SG Each set contains: scalpel, 18G needle, guidewire, dilator, catheter, 2 syringes, spigot.								
Catheter design	Catheter material	Dilator length	Dilator diameter	CH	Catheter length	Units/box	Reference	
Integrated balloon straight tip		15 cm	5.1 mm	1 0	44 cm	5	4450410	
	Cilianna	15 cm	5.9 mm	12	44 cm	5	4450412	
	Silicone	15 cm	6.6 mm	1 4	44 cm	5	4450414	
		15 cm	6.6 mm	1 6	44 cm	5	4450416	

URINARY CATHETERIZATION

Indwelling and Prostatic Catheters



NEPHROFIX®

Complete sets

	GUIDE					PIGTAIL CA	THETER		BALLOON CATHETER	
Product	Style	Tip	Ø (inch)	Ø (mm)	Length (cm)	CH8	CH11	CH14	CH12	Units/box
Nephrofix®	Seldinger	Straight	0,035"	0,9	80	4435087	4435117	4435141	-	5
Nephrofix® Certo	Schüller	J3	0,032"	0,8	90	4435583	4435613	4435648	4435745	5

All sets contain: one catheter, 3 pieces needle, guide wire, dilators, scalpel, syringe, ruler.

Exchange sets

	GUIDE					PIGTAIL CATHETE	R	BALLOON CATHETER	
Product	Style	Tip	Ø (inch)	Ø (mm)	Length (cm)	CH8	CH11	CH12	Units/box
Nephrofix® Nephrofix® Certo	Seldinger	Straight	0,035"	0,9	80	4436083	4436113	4436741	5

All exchange sets contain: one catheter, guide wire & dilator of the appropriate size.

Puncture needles: • Nephrofix®: 3 pieces 17G x 220 mm • Nephrofix® Certo: 3 pieces 18G x 220 mm - Schüller guide wire: rigid part 300 mm

H

URECATH®

Stents

Stent Tip	CLI	Length (mm)							
Stent lip	СН	200	250	275	300	*			
OPEN / CLOSE STENT DESIGN									
	5	4430921	4430930	4430948	4430956	(1)			
Tioman	6	4431057	4431073	4431081	4431090	(2)			
Tieman	7,5	4431103	4431111	4431120	4431138	(2)			
	9	-	4431316	-	-	(2)			
Round	7,5	4431200	4431219	4431227	-	(2)			

CII	Length (mm)						
CH	250			275	*		
OPEN / OPEN STENT DESIGN							
6	443749	7	4	437500	(3)		
7,5	4437535		4	437543	(3)		
	'						
CLI		Leng	th (mn	n)			
СН	200	2!	50	275	*		
CLOSE /	CLOSE / CLOSE STENT DESIGN						
7,5	4431600	4431	1618	4431626	(3)		

Catheters

Tin shana	CH	Longth (am)	Number		Eyes location	Reference	*		
Tip shape	СП	Length (cm)	of eyes	Tip	Loop	Shaft	Reference		
URETERIC CATHETERS									
	4	100	4	Х			4432690	(4)	
Straight	6	100	4	Х			4432703	(4)	
	8	100	4	Х			4432711	(4)	
	6	70	8		Х	Х	4432002	(5)	
Loop	7,5	70	8		Х	Х	4432010	(5)	
	9	70	8		Х	Х	4432029	(5)	
Loop	7	70	4		Х		4432800	(5)	

URINARY CATHETERIZATION

Urinary Dry Catheters



UROLOGICAL CATHETERS

	Length (cm)	CH	Ø (mm)	Units/box	Reference
Mann Catheter					
		■ 08	2,7	100	4500200
0		1 0	3,3	100	4500201
/7		12	4,0	100	4500202
Ш	40	1 4	4,7	100	4500203
l'i		1 6	5,3	100	4500204
1		1 8	6,0	100	4500205
		2 0	6,6	100	4500206
aton Catheter					
		■ 08	2,7	100	4500207
		1 0	3,3	100	4500208
4 3		12	4,0	100	4500209
	40	1 4	4,7	100	4500210
1		1 6	5,3	100	4500211
		1 8	6,0	100	4500212
		2 0	6,6	100	4500213
aton Catheter					
		08	2,7	100	4500214
94		1 0	3,3	100	4500215
18	10	12	4,0	100	4500216
	18	1 4	4,7	100	4500217
		1 6	5,3	100	4500218
		1 8	6,0	100	4500219

Packaging: Boxes of 10.

(1) Delivered with guide wire 900 mm x 0,6 mm (0,024"), pusher, cap, clamp and patient's card.

⁽²⁾ Delivered with guide wire 900 mm x 1,1 mm (0,043"), pusher, cap, clamp and patient's card.

⁽³⁾ Delivered with guide wire 1500 mm x 0,9 mm (0,035"), 1 side curled, pusher, cap, clamp and patient's card.

⁽⁴⁾ Delivered with introducer.

 $^{^{(5)}}$ Delivered with introducer, cap and clamp.

H

URIMED® CATH

	CH	Balloon capacity (mL)	Units/box	Reference				
URIMED® CATH FOLEY N	RIMED® CATH FOLEY NELATON							
	12	10	10	4563250				
	1 4	10	10	4563251				
	1 6	10	10	4563252				
OTTO DO	1 8	10	10	4563253				
300	2 0	10	10	4563254				
	2 2	10	10	4563255				
	2 4	10	10	4563256				
JRIMED® CATH FOLEY T	TEMANN							
	12	10	10	4563257				
	1 4	10	10	4563258				
	1 6	10	10	4563259				
	1 8	10	10	4563260				
~	2 0	10	10	4563261				
JRIMED® CATH FOLEY [DUFOUR							
	1 6	40	10	4563262				
	1 8	50	10	4563263				
	2 0	60	10	4563264				
1 Stop	2 2	70	10	4563265				
All I	2 4	80	10	4563266				

URINARY CATHETERIZATION

Indwelling Catheter Maintenance

H

URO-TAINER®

IIDA	T A I	INFR®	T
UKU	- I AI	IINEK®	ıwın

	Size	Units/box	Reference						
URO-TAINER® TWIN SUBY G									
Standard	100 mL	10	FB99839						
Twin	60 mL	10	9746609						
URO-TAINER® TW	URO-TAINER® TWIN SOLUTION R								
Standard	100 mL	10	FB99841						
Twin	60 mL (Twin)	10	9746625						

URO-TAINER® Polihexanide

	Size	Units/box	Reference
URO-TAINER® POI	IHEXANIDE		
	100 mL	10	FB99965

URO-TAINER® NaCl

	Size	Units/box	Reference			
URO-TAINER® NACL 0.9%						
	100 mL	10	FB99833			
URO-TAINER® M NACL 0.9% with injection port						
	100 mL	10	FB99853			

URINARY COLLECTION & MEASUREMENT

Urine Measurement Systems



UREOFIX® 500 CLASSIC

	Tubing length (cm)	Units/box	Reference
UREOFIX® 500 CLASSIC SETS WITH CLOSED BAGS			
Ureofix® 500 Classic with exchange bag 1,5 L	120	10	4417910
oreonx 500 classic with exchange day 1,5 L	170	10	4417920
Ureofix® 500 Classic with exchange bag 3,5 L	170	10	4417960
UREOFIX® 500 CLASSIC SETS WITH BOTTOM OUTLET BAGS			
	120	10	4417930
Ureofix® 500 Classic with premounted b.o bag 2 L	170	10	4417940
	200	10	4417950
UREOFIX® 500 CLASSIC SPARE BAGS	Volume		
Ureofix® 500 Classic closed bag	1.5 L	100	4417535
Oreotix® 500 Classic closed day	3.5 L	100	4417543
Ureofix® 500 Classic Bottom Outlet Bag	2 L	100	4417551
UREOFIX® 500 CLASSIC ACCESSORIES			
Fixation strap		10x10	4417755

H

UREOFIX® COMPACT

Production description	Tubing/Lenght (cm)	Units/Box	Reference
UREOFIX® COMPACT			
Liver Fix.® Commence 2. Chamble are 2.1 hours	150	10	4417901
Ureofix® Compact - 2 Chambers 2 L bag	120	10	4417902
Fixation strap	-	10x10	4417755



URINARY COLLECTION & MEASUREMENT

Closed Urine Drainage Bags



UREOFIX® 112 PLUS

	Units/Box	Reference
UREOFIX® 112 PLUS		
	25	284315



POST OPERATIVE 4 L BAG

Production description	Units/Box	Reference
4 L POST OPERATIVE URINE BAG		
4 L bag TUR*	20	4801020
4 L bag TUR with pump*	20	4801030

* Delivered with clamp, bed sheet clip in double pouch







URIMED® SP

	Length tube (cm)	Ø Tube (mm)	Units/box	Reference
URIMED® SP 2 L				
	100	9	10	28610
-	100	9	50	28650
URIMED® SP+ 2 L				
	130	9	30	28630

URINARY COLLECTION & MEASUREMENT

Adult Urine Bags



URIMED® BAG

Night and Day bags (2 L)

	Length tube	Ø Tube (mm)	Units/box	Reference		
URIMED® BAG 2 L						
Non sterile	120 cm connected connector	9	10	28300A		
MEDICARE® 2 L						
Sterile	90 cm	6.5	10	28115		
Non sterile	90 cm	6.5	10	28114		

Leg bags

	Length tube	Ø Tube (mm)	Units/box	Reference		
URIMED® BAG F	PLUS 0,5 L					
Non sterile	60 cm cut to fit	9	10	28501A		
URIMED® TRIBAG PLUS 0,8 L						
3 Chamber – Sterile	60 cm connected connector	9	10	28305		
Non sterile	60 cm cut to fit	9	10	28306		
URIMED® BAG PLUS 1.5 L						
Non sterile	90 cm cut to fit	9	10	28150A		

URIMED® ACCESSORIES

Accessories

Tube		Units/box	Reference				
URIMED® FIX*							
	Small	2	68520R				
	Medium	2	68530R				
	Large	2	68540R				
URIMED® STRAPS							
		1 pair	68550R				
* Distributed by	B. Braun	•					

	Tube	Units/box	Reference				
URIMED® KLETT							
	Small	10	28260				
	Large	10	28261				
URIMED® HANDLE							
		1	68548R				





URIMED® B'BAGS

Closed bags

		Capacity (L)	Units/box	Reference			
URIMED® B'BAG – Tube length: 90 cm – Tube diameter intR : 5.05 & extR 6.5							
Sterile	No ARV*	2	100	29400			
bags	INO AITV	1,5	100	29420			
Non sterile bags	No ARV*	2	250	29415			
		1,5	250	29435			
	With ARV*	2	250	29410			
		1,5	250	29430			
Sterile bags	No ARV*	2	30 + 1 hanging accessory**	29416			
Non sterile bags	With ARV*	2	30 + 1 hanging accessory**	29411			

Drainable bags

		Capacity (L)	Units/box	Reference			
URIMED® B'BAG - Tube length: 90 cm - Tube diameter intR : 5.05 & extR 6.5							
Sterile bags	With ARV*	2	100	29440			
Non sterile	With ARV*	2	200	29450			
bags	VVILII ANV	1,5	200	29455			
Non sterile bags	With ARV*	2	30 + 1 hanging accessory**	29451			

^{*} ARV= Anti- Reflux valve **made of PVC: DEHP free

URINARY COLLECTION & MEASUREMENT

Pediatric Urine Bags





URINOCOL®

		Units/box	Reference
URINOCOL PAEDIATRIC			
Sterile ope	en Boy	50	227546*
paediatric	bags Girl	50	227566*
Sterile clos	sed Boy	100	227550*
paediatric	bags Girl	100	227560*

^{*}For each country, a specific letter is added at the end of the reference.

H

STANDARD PEDIATRIC URINE BAGS To be confirmed

	Sterile/ Non sterile	ARV */ Without ARV*	Closed / Drainable	Units/box	Reference
STANDARD PEDIA	TRIC URINE BAGS				
	Non sterile	Without ARV*	Drainable - With drainage tube of 90 cm - With slide clamp	20 – Individually packed	4423046
For Children	Non sterile	With ARV*	Drainable - With bottom outlflow nozzle	10 - Bulk packed	4423135
	Sterile – Individually wrapped	Without ARV*	Closed	10 - Individually packed	4425030
	Non sterile	Without ARV*	Closed	10 - Bulk packed	4423038
For Neonates	Sterile – Individually wrapped	Without ARV*	Closed	10 - Individually packed	4425049
ARV= Anti-reflux va	ve				

Weight (g) Units/box Reference Standard 1200-2500 50 227571* Sterile premature bags Mini < 1200 227570* 50 Sterile prema-Standard 1200-2500 50 227573* ture bags (with extension line**) 50 227572* Mini <1200

^{**}made of PVC: DEHP free

BOWEL MANAGEMENT & URINARY INCONTINENC



RECTOBAG®

	Product de	escription	Units/box	Reference
RECTOBAG® FLUID COLLECTING	g systems			
	Closed System for intestinal I	avage		
	2 L lavage bag	Rectal tube with balloon		
	150 cm large bore tubing with side clamp	Omnifix® 50/60 mL Luer Lock syringe	20	4423364
1 Name 1	Y-piece with connector and side clamp			
	Sterile / single-packed			
•	Closed System for intestinal lavage and irrigation			
illo	• 2 L lavage bag	• 150 cm large bore tubing with side	20	4423348
T	Connector	clamp	20	
0	Sterile / single-packed			
	Closed System for lavage bag			
			40	4423356
18/11-20	Non sterile / bulk packed			

BOWEL MANAGEMENT & LIRINARY INCONTINENCI

Rectal Catheters



RECTAL CATHETERS

	Length (cm)	CH	Units/box	Reference
RECTAL CATHETERS				
		2 2	100	4500300
	40	■ 25	100	4500301
	40	2 8	100	4500302
		■ 30	100	4500303

Ø Tube (mm) Units/box Reference

BOWEL MANAGEMENT & URINARY INCONTINENC

Male External Catheters



URIMED® VISION

URIMED® VISION STANDARD – 7,5 cm length with 4,5 cm adhesive area					
	25	30	IH2525		
	29	30	IH2529		
	32	30	IH2532		
	36	30	IH2536		
	41	30	IH2541		
URIMED® VISION ULTRA	URIMED® VISION ULTRA – 7,5 cm length with full adhesive area				
11	25	30	IH3525		
	29	30	IH3529		
	32	30	IH3532		
	36	30	IH3536		
		30	IH3541		

	Ø Tube (mm)	Units/box	Reference		
URIMED® VISION SHORT – 4,5 cm length with full adhesive area					
m .	25	30	IH4525		
	29	30	IH4529		
	32	30	IH4532		
	36	30	IH4536		
	41	30	IH4541		

	Description	Units/box	Reference
URINE BAGS			
	Urimed® Bag Plus 500 mL	30	28501A
	Urimed® Bag Plus 800 mL	10	28306A
	Urimed® Bag Plus 1500 mL	30	28150A
	Urimed® Bag 2000 mL	30	28300A

BOWEL MANAGEMENT & URINARY INCONTINENCE

Female Intravaginal Tampon



DIVEEN®

DIVEEN®		Size	Units/box	Reference
-58		Small	5 Devices + 1 Applicator	3132
	1000	Medium	5 Devices + 1 Applicator	3133
and the same	8	Small	15 Devices + 1 Applicator	3022
		Medium	15 Devices + 1 Applicator	3023

This document, its contents, including institutional data, information, trademarks and logos mentioned here are the exclusive property of B. Braun. Any representation and / or reproduction, total or partial, of this document and its contents without the express prior consent of B. Braun, is strictly prohibited and constitutes an infringement of the intellectual property rights of B. Braun.

Commercialized by:

B. Braun Medical I 204 avenue du Maréchal Juin I 92100 Boulogne Billancourt - France Tél. 01 41 10 53 00 l Fax 01 41 10 53 99 l www.bbraun.fr

B. Braun Medical SAS with a share capital of 31 000 000 € I RCS Nanterre 562050856



For most countries a specific letter is added at the end of the radical reference.

All volumes are measured according to ISO 8670-2 standard.
This medical device is a regulated health care product and on the basis of this regulation is certified CE.

Non-biding documents and photographs.

The products commercialized by B. BRAUN MEDICAL shall exclusively be used in accordance with the package inserts which have been updated and included in the