



2-piece system

FLEXIMA[®] KEY

THE KEY TO FEELING SECURE WITH EXTRA FLEXIBILITY

Flexima[®] Key

The guided positioning

A major step towards security

Key guided positioning

- Precise and simple positioning with a guiding system
- Minimum patient education: intuitive coupling



Key guiding system

Secured coupling system


- Large stable track, no distortion.
- Specific adhesive foam formula



Bag fitting

Easy bag removal

- Special cut-away section at the top of the flange
- Bag ring leaves no residue when the bag is removed

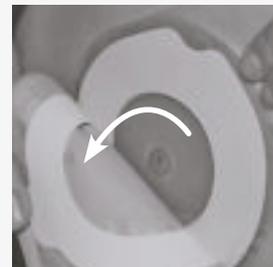

Removal

Flexibility

Soft material for flexibility

- Follows body contours

Smooth and low profile track


Soft convexity

- For flush or slightly retracted stomas
- A belt can be used




Flexima[®] Key Drainable Roll'Up bags

High performance filter

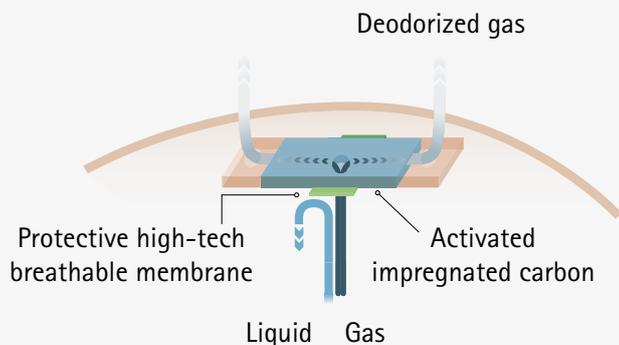

- Leak-Proof and permeable to gases high-tech breathable membrane

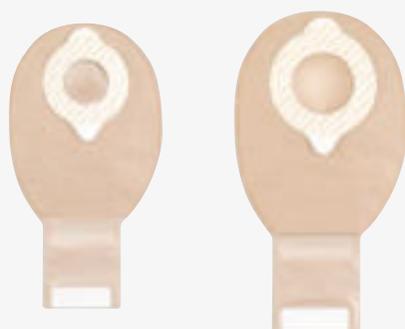

Midi:
525 ml*




Maxi:
610 ml**

Shape and features

- Ergonomic shape
- Optimized sizing

Integrated soft Roll'Up clamp

Intuitive and secure closing system

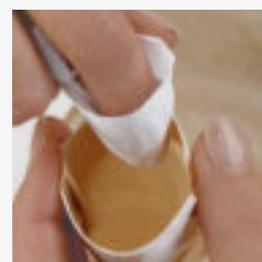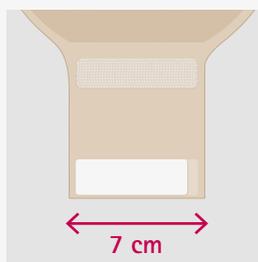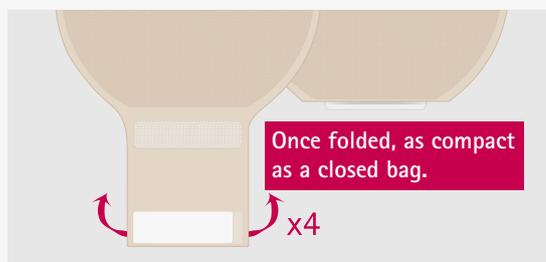

- Soft edges

- A large 7 cm outlet for both colostomates and ileostomates


- White opening tabs

* for diameter 50mm
**for diameter 80mm

Flexima[®] Key

Closed, drainable and uro bags

Closed bag

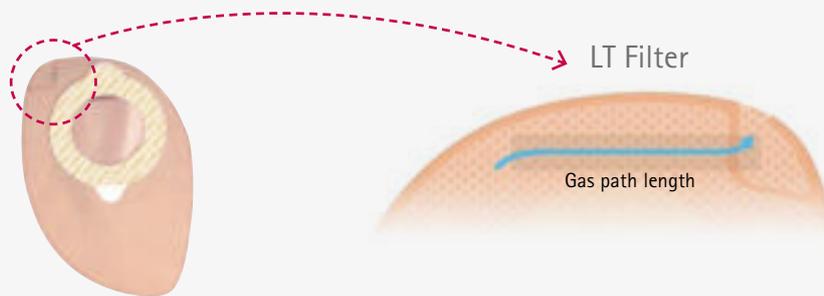


- **Laminar technology filter:**
- Optimized performance of degassing and deodorizing

Maxi drainable bag



- Funnel-shaped
- Large diameter direct connection to a Flow Collector bag
- Suitable for patients with low dexterity

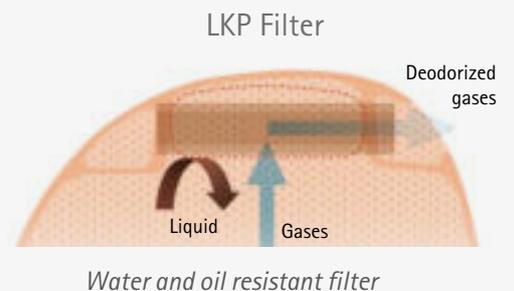

Flexima[®] Key Uro



A Secure outlet OUTLET FOR DAY

- Position open/closed
- Fits in perfectly strengthened by a click and a final locking system
- Easy handle; no risk of wet fingers

OUTLET FOR NIGHT

- Click at connection
- Can be connectd to a night drainage bag
- Outlet designed to avoid bag twist



Closed



Open



Universal Connection



Capacity: 2L

Flexima[®] Key

Right positioning – Fit from the start



1 Insert the tab into the wafer slot

2 Stick gently from bottom to top

Users very satisfied*

- Security feeling
- Comfort during wear

* Report on Flexima[®] Key post-marketing user evaluation - 2005

« The key system is «idiot proof; always in the right spot»
Giny- ET nurse, Netherlands

« Very easy to use for a lot of people»
Jos- ET nurse, Netherlands

« The Key is a good idea because the problem we're facing with this kind of product is being able to position it correctly.
With Key, you can't go wrong!»
Isabelle- ET nurse, Netherlands

Flexima® Key

FLEXIMA® KEY FLAT

✂: BOX OF 5 / PRE-CUT: BOX OF 5

	∅ (mm)	References
✂	40/ 15-35	62040
✂	50/15-45	62050
✂	60/ 15-55	62060
✂	80/ 15-75	62080

FLEXIMA® KEY CONVEX

✂: BOX OF 5

✂	50/ 15-25	62052
✂	60/ 15-35	62062
✂	60/ 15-45	62063
	60/ 25	62064
	60/ 30	62066
	60/ 35	62067
	60/ 40	62069

FLEXIMA® KEY CLOSED - BOX OF 30

Size	∅ (mm)	References	Volume
Midi	40	● 62140	±440ml
		○ 62240	
	50	● 62150	±420ml
		○ 62250	
	60	● 62160	±550ml
		○ 62260	
	80	● 62180	±520ml
		○ 62280	

✂ Cut-to-fit | ● Beige | ○ Transparent

For most countries a specific letter is added at the end of the radical reference

FLEXIMA® KEY ROLL'UP - BOX OF 30

Size	∅ (mm)	References	Volume
Midi	40	● 621241	±550ml
		○ 621240	
	50	● 621251	±525ml
		○ 621250	
Maxi	60	● 621261	±670ml
		○ 621260	
	80	● 621281	±610ml
		○ 621280	

FLEXIMA® KEY URO - BOX OF 30*

Midi	40	● 62740	±580ml
		○ 62840	
	50	● 62750	±570ml
		○ 62850	
	60	● 62760	±570ml
		○ 62860	

*with 30 universal connectors

FLEXIMA® KEY MAXI DRAINABLE- BOX OF 30

Maxi	50	● 62550	±780ml
		○ 62650	
	60	● 62560	±730ml
		○ 62660	
	80	○ 62680	±655ml

FLOW COLLECTOR®

+	039900	BOX OF 15	±2l
	039901	BOX OF 5	±2l

Stomacare belt	References
1 unit	5028

This document, its contents, including institutional data, information, trademarks and logos mentioned herein are the exclusive property of B. Braun. Any representation and / or reproduction, total or partial, of this document and its contents without the express prior consent of B. Braun, is strictly prohibited and constitutes an infringement of the intellectual property rights of B. Braun.

Commercialized and manufactured by :

B. Braun Medical SAS | 26 rue Armengaud | 92210 Saint-Cloud - France

Tél. 01 41 10 53 00 | Fax 01 70 83 45 00 | www.bbraun.fr | www.bbraun.com

Société par Actions Simplifiée with a share capital of 31 000 000 € | RCS Nanterre 562050856



All stoma bags volumes are measured according to ISO 8670-2 standard.
This medical device is a regulated health care product and on the basis of this regulation is certified CE.
The products commercialized by B. BRAUN MEDICAL shall exclusively be used in accordance with the package inserts which have been updated and included in the boxes.
Non binding documents and photographs

www.bbraun.com/living-with-a-stoma