

# Omnitest<sup>®</sup> 5

## BLOOD GLUCOSE MONITORING SYSTEM

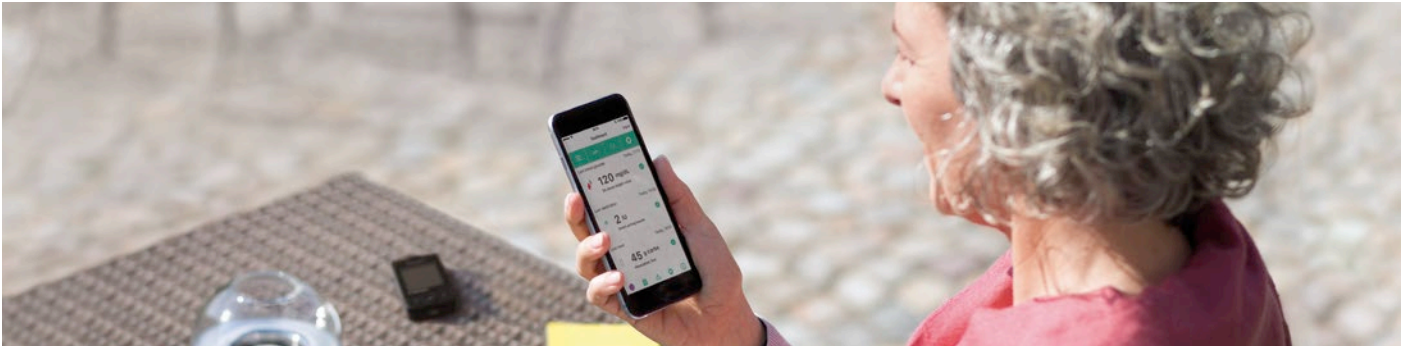


Image in original size

Omnitest<sup>®</sup> 5 blood glucose meters convince with a combination of modern design, easy handling and accurate measurement results.


- 0.5 µL sample volume
- 5 seconds measuring time
- Insufficient blood sample alert
- Marking function
- 5 standard alarms and timer
- Average results for 7, 30 and 90 days
- Memory capacity of 500 readings
- Electronic self-test prior each measurement
- Automatic data transfer via Bluetooth<sup>®</sup> technology to Omnitest<sup>®</sup> app or via USB to DIABASS<sup>®</sup>, Accu-Chek<sup>®</sup> Smart Pix or diasend<sup>®</sup>. More info: [www.bbraun.com/diabetes-datamanagement](http://www.bbraun.com/diabetes-datamanagement)

### MEETS THE REQUIREMENTS OF EN ISO 15197:2015




[www.bbraun.com/omnitest5](http://www.bbraun.com/omnitest5)

Measuring range	10 – 600 mg/dL
Hematocrit	20 – 60%
Enzyme	FAD glucose dehydrogenase
Blood sample	Capillary or venous whole blood
Calibration	Plasma
Code	Automatic code recognition
Dimensions (L x W x H)	78 mm x 43 mm x 16 mm
Power supply	Two 3V lithium batteries (CR2032)

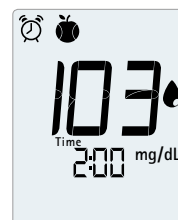

#### TIMER

Beside 5 standard alarms with a fixed time, there is also a timer function: The user can set a reminder of a follow-up measurement (e.g. 30, 60 or 120 min. after a meal) just at the push of a button.

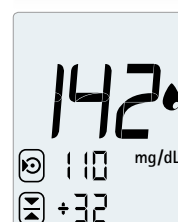

#### TARGET LEVEL FUNCTION

Omnitest<sup>®</sup> 5 does not only show the measurement result but also its deviation from the pre-defined target level. This helps to better understand the result.

# Omnican<sup>®</sup> Lance

## LANCING DEVICE AND LANCETS



[www.bbraun.com/omnican-lance](http://www.bbraun.com/omnican-lance)

Images in original size

Thanks to its easy handling, the new lancing device Omnican<sup>®</sup> Lance offers more comfort during blood collection.

- Less painful due to straight and fast lancet leading
- 7 finely graduated depth settings for different skin types
- One-hand operation
- Coloured signal upon loading

In combination with Omnican<sup>®</sup> Lance soft lancets, the new lancing device is your reliable partner for a soft blood collection.

- Gentle puncture and ideal blood flow
  - silicone oil coated tip
  - three-facet grinding
- Maximum comfort thanks to thin 30G (ø 0.30 mm) or 28G (ø 0.36 mm) lancet tip

